Body Armor[®] Heel Reliever


- EN** Heel Protector
- DE** Fersen-Freilagerung
- ES** Soporte abierto para el talón
- FR** Dispositif de décharge du talon
- IT** Supporto aperto per tallone
- NL** Hielbeschermer
- DA** Hælsuspension
- SV** Hälavlastare
- PL** But odciążający piętę
- SL** Opornica za sprostitev pete
- HR** Ortoza za rasterećenje pete
- PT** Descanso do calcanhar





DARCO

Body Armor® Heel Reliever

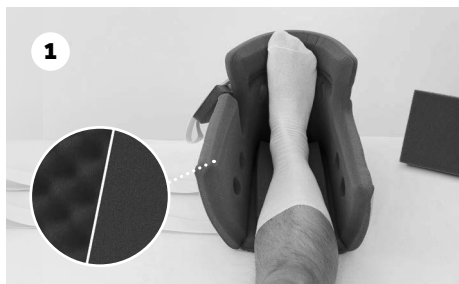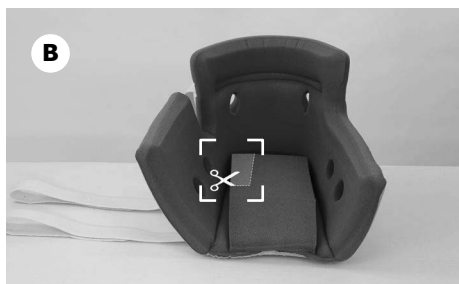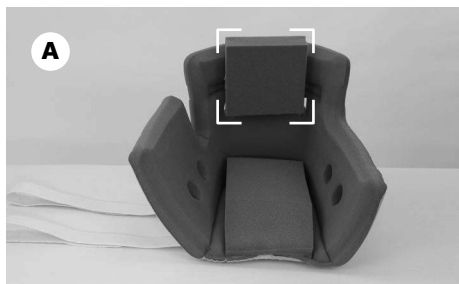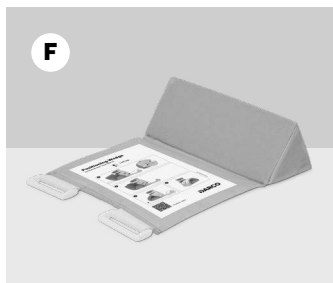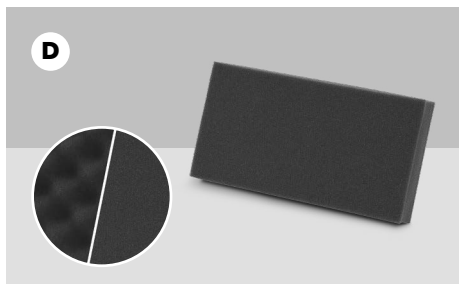Body Armor® Heel Reliever | Heel Protector

Important notes

If there is a serious deterioration in your state of health as a result of using the medical device, please report this to your doctor, specialist dealer or the manufacturer as soon as possible. If you report this to the specialist dealer, the dealer will forward this immediately. The reports are agreed with and processed by the competent authorities.

Please always follow these instructions for use. In case of adverse effects (e. g. skin swelling, pressure points, disorders of blood circulation etc.) please tell your doctor. Function and optimum protection are no longer guaranteed if the product is handled incorrectly. This aid assists recovery and should be used exclusively under the instruction of and in consultation with your doctor or a healthcare professional. DARCO (Europe) GmbH does not assume any responsibility. Do not wear the product with injured, irritated, damaged or intact skin. You should use bandages or socks, for example.

Note: This device is allowed for single-patient use only. It does not need to be cleaned or reprocessed between uses. The product is only intended for support. Patients should not walk when wearing the product.

Purpose

The Body Armor® Heel Reliever heel protector is for the prevention and healing of heel pressure sores as well as diabetic foot syndrome and decubitus.

Indications / Contraindications

Indications

- For the prevention and healing of heel pressure ulcers as well as diabetic foot syndrome and decubitus
- Braden scale scores of 14 or less
- Diabetes with neurological sensory deficit
- Pedal pulse not palpable
- Paraplegia
- Prolonged immobility or coma

Contraindications

- Wounds on the calf
- Spasticity
- Contraction of the ankle

Design Features

The Body Armor® Heel Reliever heel protector is available in 2 versions: Convoluted or smooth. Both product variants are laminated with a smooth fabric cover to allow the leg to „glide“ without changing the position of the exposed area. The correct fixation of the product on the patient is achieved with three straps that hold the forefoot and the lower leg in an orthograde position. If the leg is rotated excessively over the hip joint, the desired position of the leg can be achieved with the positioning wedge. An extra pad for the calf area and dorsal foot area (included in the package) enables individualised pressure redistribution and safe application of the product.

Use

1. Open the Body Armor® Heel Reliever heel protector. Remove the extra foam pad (photo D). Set the extra pad aside and keep it available. You may need it later to customise the fit.
2. Place the patient's foot inside the support with the heel positioned above the heel suspension opening. The heel should be positioned below the fixed foam pad (photo 1 + 2).
3. Pull the hook and loop straps across the shin towards the loops. Guide the straps through the loops. Secure the straps with the hook and loop fasteners (photo 3). Allow the side flaps to overlap if the calves are particularly narrow.
4. Check the fit. There should be sufficient space for your fingers between the heel opening and the bed (photo 4). Do not fasten the heel protector too loosely or too tightly around the leg. The hook and loop straps should not be in contact with the skin. The forefoot pad (photo E) can be used optionally for protection between the hook and loop strap and foot (photo 5).
5. Instructions for use: Depending on the patient / indication, either use a heel protector with smooth or convoluted surface (photo 1).
6. Instructions for use: The Body Armor® Heel Reliever can be used without restriction in beds with adjustable lying surfaces as well as in beds with side rails.

Customising the fit

- Foot drop: Secure the extra foam pad in a vertical position behind the sole of the foot (photo A).
- Malleolar decubitus (ankle ulcer): Cut away the nubs surrounding the ankle on the inside of the support. Cut a portion of the fixed foam pad if necessary (photo B).

- Hip rotation: Fix the positioning wedge (photo F) to the heel protector, according to use.
- Bariatric legs: If the heel comes into contact with the bed sheet, add the extra pad on top of the fixed pad for increased elevation (photo C).
- Achilles tendon: If redness occurs in this area, cut a "V" shape out of the fixed pad. If using both foam pads, also customize the upper foam pad (photo C).

Note: The pressure on the heel is redistributed to the heel cord and gastrosoleus complex. The product should be removed regularly according to local guidelines and the leg and foot should be examined for skin lesions. Should redness (erythema) occur, a V-shaped wedge should be removed from the fixed pad in the area of the heel cord. Oedematous patients should use the smooth Body Armor® Heel Reliever version to avoid skin irritation.

For more information on the current EPUAP guidelines, please scan:

Combination with other products

No specifications

Working life / Lifespan of the product

The lifespan of the medical device is determined by natural wear and tear due to proper handling / cleaning. The medical device is designed to be used for the usual duration of use prescribed by a doctor without showing a reduction in safety or performance. After the period of use prescribed by your doctor, check the function of the product together with your doctor or specialist dealer.

Technical data

Size

Universal (25 – 54 cm cm calf circumference, up to shoe size 47)

Universal fit for right and left

Colour

Grey-Blue

Material

Heel protector

Outer cover: Polyester

Foam: Latex-free polyurethane foam

Hook and loop straps: Nylon

Forefoot pad: Polyester, latex-free polyurethane foam

Extra foam pad: Latex-free polyurethane foam

Positioning wedge

Cover: Polyester

Wedge: Ethylene vinyl acetate copolymer

Laundry bag: Polyester

Weight

Convuluted: 512 g

Smooth: 496 g

Maintenance

The product does not require maintenance.

Care instructions

Hand wash

- Use mild detergent and warm water. Allow to dry at room temperature.

Machine wash

- Close hook and loop straps and wash in the provided laundry bag (photo G) or a pillowcase (60° C). Allow to dry at room temperature.
- Disinfection after washing: With an antiseptic spray

Water, standard cleaning agents and disinfectants can be used for cleaning. The hook and loop strips should always remain closed during cleaning.

Disposal

Dispose of after use. Please dispose of the packaging in an environmentally friendly manner. Dispose of this product after use in accordance with applicable local laws and regulations.

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ALBO-Healthcare GmbH
DARCO UK Ltd